

Calcium Chloride Solution

CaCl₂ SOLUTION

I Revision bar indicates update to previous version.



Intended Use

Calcium Chloride Solution is used as Supplementary Reagent for various coagulation assays.

Reagents

Reagent	Description	Storage	Stability
Calcium Chloride Solution CaCl ₂ SOLUTION	Ready to use liquid containing: • Calcium chloride (0.025 mol/L)	2–25 °C May be used up to the expiry date indicated on the label if stored unopened.	2–25 °C: once opened, 8 weeks ^a

^a closed original vial

On-board stability

Information regarding on-board stability is specified in the Reference Guides (Application Sheets) for the different coagulation analyzers.

Warnings and Precautions

For *in-vitro* diagnostic use only.

For laboratory professional use.

According to EU regulation 2017/746, any serious incident that has occurred in relation to the device shall be reported to the manufacturer and the competent authority of the EU Member State through your local distribution representative in which the user and/or patient is established.

Safety data sheets (MSDS/SDS) available upon request.

Dispose of hazardous or biologically contaminated materials according to the practices of your institution. Discard all materials in a safe and acceptable manner and in compliance with all government requirements.

Procedure

Materials Provided

REF	Contents		
ORHO37	Calcium Chloride Solution CaCl ₂ SOLUTION	10 ×	15 mL

- CaCl_2 **SOLUTION** is added automatically by the respective coagulation analyzer when performing the respective assay.
- Please refer to the respective Instructions for Use and Application Sheets for the assays to be used with CaCl_2 **SOLUTION**.

Materials Required but not Provided

Item	Description
-	Reagents and associated materials according to the respective Instructions for Use of the coagulation assays to be used with CaCl_2 SOLUTION .
Coagulation analyzers ^b , such as:	• Automated Blood Coagulation Analyzer CA-600 series (CA-600 series) • AUTOMATED BLOOD COAGULATION ANALYZER CS-2500 (CS-2500 System) • AUTOMATED BLOOD COAGULATION ANALYZER CS-5100 (CS-5100 System) • Automated Blood Coagulation Analyzer CN-3000/CN-6000 (CN-3000/CN-6000 System)

^b Availability of analyzers may vary by country.

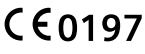
Please note that the applications on other analyzers can be validated by the instrument manufacturer in accordance with the requirements of the REGULATION (EU) 2017/746 under their responsibility as long as the intended purpose and performance are not modified.

Technical Assistance

For customer support, contact your local technical support provider or distributor.

Definition of Symbols

The following symbols may appear on the product labeling:

	Do not reuse		Use By
	Batch Code		Catalogue Number
	Caution		Manufacturer
	Authorized representative in the European Community		Authorized representative in Switzerland
	Contains sufficient for <n> tests		Biological Risks
	<i>In Vitro</i> Diagnostic Medical Device		Temperature Limitation
	Consult instruction for Use		Non-sterile
	CE marking of conformity		CE marking of conformity with notified body ID number. Notified body ID number can vary.
	Contents		Reconstitution volume
	Level		Keep away from sunlight and heat
	Warning		Danger
	Prescription device (US only)		Device Identification (UDI) barcode
	REACH Authorization Number xx/xx/xx		

Legal Information

Sysmex is a trademark of SYSMEX CORPORATION.

All other trademarks are the property of their respective owners.

© Siemens Healthineers, 2010–2024. All rights reserved.



Siemens Healthcare Diagnostics Products GmbH

Emil-von-Behring-Str. 76

35041 Marburg

Germany

[siemens-healthineers.com](https://www.siemens-healthineers.com)